

PIPETTE TIPS

LABSOLUTE® offers a complete assortment for your pipetting tasks – no matter whether packed in bags, within racks or as an efficient refill system: all pipette tips are part of a system in which the components can be exchanged. The pipette tips are made of pure PP and comply with EN ISO 8655. In combination with the LABSOLUTE® pipettes, they constitute the perfect pipetting unit.

Thanks to their universal fit and the high-quality material, the tips guarantee a firm and absolutely tight hold on various pipette models. And all this with much less attachment and ejection forces. The fully automated production and packaging of the pipette tips under clean room conditions guarantee completely contamination-free working conditions with products which are certified free of DNase, RNase, human DNA and pyrogens – to ensure reproducible and consistent results.

- High quality standard
- Graduated tips: the volume is visible at first glance
- 200 and 1000 µl tips are available in yellow or blue, respectively
- Autoclavable

Our pipette tips boxes for comfortable and efficient pipetting:


- Simple-to-use closure for single-handed opening
- The transparent design provides a complete overview of the tips
- Coloured trays enable fast recognition of the volume
- Refill system: easy and fast removal and replacement of the trays
- Stackable boxes for space-saving and safe storage
- Autoclavable at 121 °C



PIPETTE TIPS WITHOUT AEROSOL BARRIER IN BAGS, NON-STERILE

The 0.1–10 XL / 20 µl tips (7.695 841) are compatible with 20 µl Eppendorf Reference®, Eppendorf Research® Plus and Brand Transferpette®. Bevelled tips (7.695 842 and 7.695 843): due to the bevelled end the contact surface of the pipetted liquid is reduced and less residual liquid remains at the tip. Item 7.695 856 is delivered in a useful folded box.

Volume range µl	Figure	Colour	Overall length mm	Grad. µl	PK	Art. no.
0,1–10	1	Transparent	31	2, 10	1000	7.695 840
0,1–10 XL/20	2	Transparent	44	1, 5, 10, 20	1000	7.695 841
2–200 (bevelled tip)	3	Transparent	51	-	1000	7.695 842
2–200 (bevelled tip)	4	Yellow	51	-	1000	7.695 843
2–200	5	Transparent	53	10, 50, 100	1000	7.695 844
2–200	6	Yellow	53	10, 50, 100	1000	7.695 845
20–300	7	Transparent	55	50, 100, 200	1000	7.695 846
100–1000	8	Transparent	71	100, 500	500	7.695 847
100–1000	9	Blue	71	100, 500	500	7.695 848
100–1250	10	Transparent	76	-	500	7.695 849
100–1250 XL	11	Transparent	98	200, 500, 1000	500	7.695 854
500–5000, Eppendorf Fit	12	Transparent	120	1000, 1500, 2000, 2500, 3000, 3500, 4000, 4500, 5000	250	7.695 850
500–5000, Gilson Fit	13	Transparent	120	1000, 1500, 2000, 2500, 3000, 3500, 4000, 4500, 5000	250	7.695 933
500–5000, Thermo Fit	14	Transparent	147	-	250	7.695 855
500–5000, Thermo Fit	14	Transparent	147	-	100	7.695 856
1000–10000, Eppendorf/ Gilson Fit	15	Transparent	160	1000, 1500, 2000, 2500, 3000, 3500, 4000, 4500, 5000, 5500, 6000, 6500, 7000, 7500, 8000, 8500, 9000, 9500, 10000	200	7.695 934
1000–10000, Eppendorf Fit	16	Transparent	165	1000, 5000, 10000	100	7.695 851
2000–10000, Thermo Fit	17	Transparent	150	-	100	7.695 857

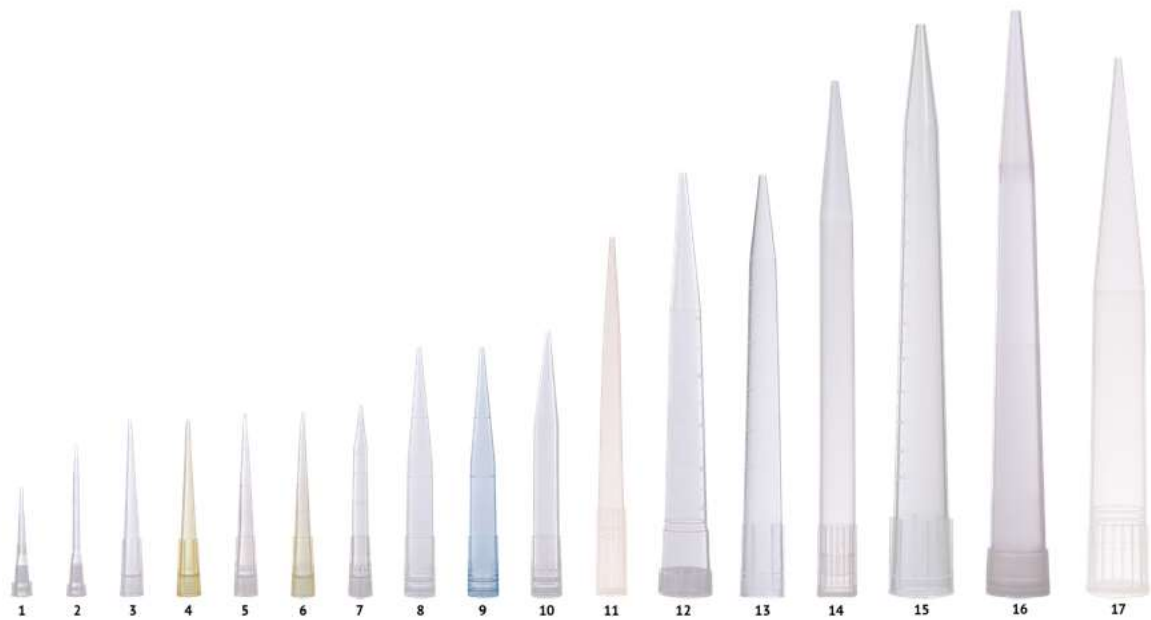

PIPETTE TIPS WITHOUT AEROSOL BARRIER IN RACKS, NON-STERILE

The 0.1–10 XL/20 µl tips (7.695 861) are compatible with 20 µl Eppendorf Reference®, Eppendorf Research® Plus and Brand Transferpette®.

Volume range µl	Figure	Colour	Overall length mm	Grad. µl	Pack size pieces	PK	Art. no.
0.1–10	1	Transparent	31	2, 10	10 x 96	960	7.695 860
0.1–10 XL/20	2	Transparent	44	1, 5, 10, 20	10 x 96	960	7.695 861
2–200	5	Transparent	53	10, 50, 100	10 x 96	960	7.695 862
2–200	6	Yellow	53	10, 50, 100	10 x 96	960	7.695 863
20–300	7	Transparent	55	50, 100, 200	10 x 96	960	7.695 864
100–1000	8	Transparent	71	100, 500	10 x 96	960	7.695 865
100–1000	9	Blue	71	100, 500	10 x 96	960	7.695 866
100–1250	10	Transparent	76	-	10 x 96	960	7.695 867
100–1250 XL	11	Transparent	98	200, 500, 1000	10 x 96	960	7.695 868

PIPETTE TIPS WITHOUT AEROSOL BARRIER IN RACKS, STERILE

The 0.1–10 XL/20 µl tips (7.695 871) are compatible with 20 µl Eppendorf Reference®, Eppendorf Research® Plus and Brand Transferpette®.


Volume range µl	Figure	Colour	Overall length mm	Grad. µl	Pack size pieces	PK	Art. no.
0.1–10	1	Transparent	31	2, 10	10 x 96	960	7.695 870
0.1–10 XL/20	2	Transparent	44	1, 5, 10, 20	10 x 96	960	7.695 871
2–200	5	Transparent	53	10, 50, 100	10 x 96	960	7.695 872
2–200	6	Yellow	53	10, 50, 100	10 x 96	960	7.695 873
20–300	7	Transparent	55	50, 100, 200	10 x 96	960	7.695 874
100–1000	8	Transparent	71	100, 500	10 x 96	960	7.695 875
100–1000	9	Blue	71	100, 500	10 x 96	960	7.695 876
100–1250	10	Transparent	76	-	10 x 96	960	7.695 877
100–1250 XL	11	Transparent	98	200, 500, 1000	10 x 96	960	7.695 878

FILTER TIPS IN RACKS, STERILE

The LABSOLUTE® filter tips offer dual protection: the fine-pored filter made of PE serves as a defence to protect the samples from aerosols and biomolecules. At the same time, the hydrophobic material prevents the pipette from contamination.

The 0.1–10 XL/20 µl tips (7.695 881) are compatible with 20 µl Eppendorf Reference®, Eppendorf Research® Plus and Brand Transferpette®.

Volume range µl	Figure	Colour	Overall length mm	Grad. µl	Pack size pieces	PK	Art. no.
0.1–10	17	Transparent	31	2, 10	10 x 96	960	7.695 880
0.1–10 XL/20	18	Transparent	44	1, 5, 10, 20	10 x 96	960	7.695 881
0.5–20	19	Transparent	53	20	10 x 96	960	7.695 882
2–100	20	Transparent	55	10, 50, 100	10 x 96	960	7.695 883
2–200	21	Transparent	55	10, 50, 100, 200	10 x 96	960	7.695 884
20–300	22	Transparent	55	50, 100, 200	10 x 96	960	7.695 885
100–1000	23	Transparent	76	-	10 x 96	960	7.695 886
100–1250 XL	24	Transparent	98	200, 500, 1000	10 x 96	960	7.695 887

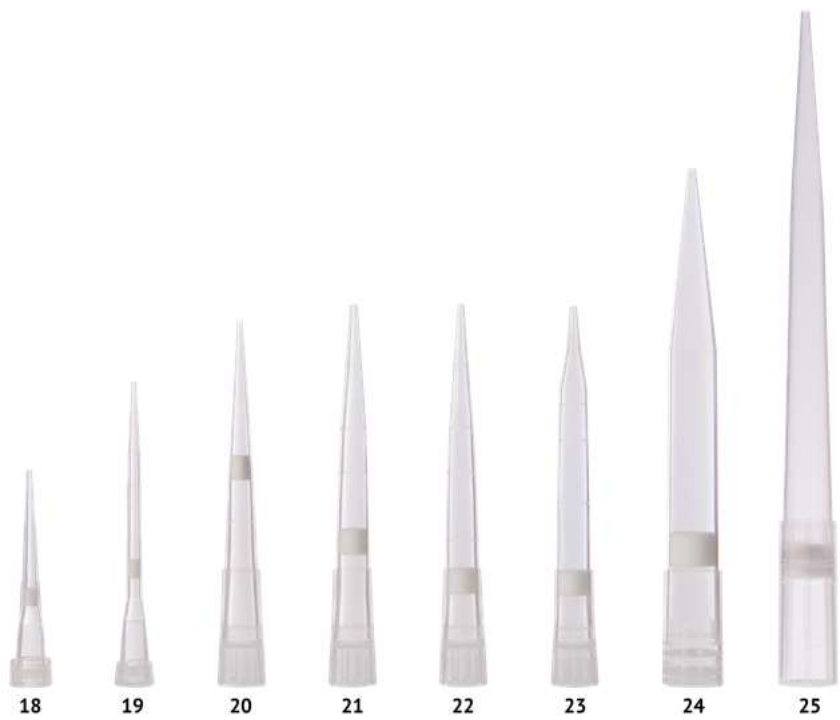

REFILL SYSTEMS, PIPETTE TIPS WITHOUT AEROSOL BARRIER, NON-STERILE

The environmental and budget friendly LABSOLUTE® pipette tip boxes provide trays which can be easily replaced and directly indicate the respective volume due to their colour. To safely store the pipette tips, the refill systems are delivered in the form of two individually shrink-wrapped towers. The outer packaging made of cardboard additionally guarantees that the pipette tips remain dust-free. Last but not least each refill system contains **two refillable pipette tip boxes for free!**

The 0.1–10 XL/20 µl tips (7.695 891) are compatible with 20 µl Eppendorf Reference®, Eppendorf Research® Plus and Brand Transferpette®.

Volume range µl	Figure	Colour	Overall length mm	Grad. µl	Pack size pieces	PK	Art. no.
0.1–10	1	Transparent	31	2, 10	10 x 96	960	7.695 890
0.1–10 XL/20	2	Transparent	44	1, 5, 10, 20	10 x 96	960	7.695 891
2–200	5	Transparent	53	10, 50, 100	10 x 96	960	7.695 892
20–300	7	Transparent	55	50, 100, 200	10 x 96	960	7.695 893
100–1000	8	Transparent	71	100, 500	10 x 96	960	7.695 894

PIPETTE TIPS BOXES

The unfilled trays of the LABSOLUTE® refillable pipette tips boxes hold up to 96 tips and are coloured differently according to the volume of the appropriate tips. The trays are autoclavable at 121 °C and can be replaced easily as it is the case for all LABSOLUTE® pipette tips boxes. The boxes are suitable for the LABSOLUTE® pipette tips (7.695 840–7.695 894).

Description	Suitable for	PK	Art. no.
Pipette tips boxes	10 µl, 10 µl XL/20 µl	10	7.695 897
Pipette tips boxes	20 µl–300 µl	10	7.695 898
Pipette tips boxes	1000 µl, 1250 µl	10	7.695 899

RACKS FOR 5 ML AND 10 ML PIPETTE TIPS

The LABSOLUTE® boxes made of PP offer space for 18 pieces of 5 ml (7.695 850) or 10 ml (7.695 851) pipette tips from LABSOLUTE®, respectively. Tips are not included.

Suitable for	Height Box mm	Ø mm	PK	Art. no.
5 ml pipette tips	188	102	2	7.695 852
10 ml pipette tips	188	102	2	7.695 853

